



# BioGX

Molecular Made Easy

## Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

**REF** 450-008-Series



24 Reactions 450-008-C-MAX




64 Reactions 450-008-LMP

## Product Insert

For Research Use Only: Not intended for In Vitro Diagnostic Use

For use with BD MAX™ System, ABI QuantStudio™ 5, Bio-Rad CFX96  
Touch™

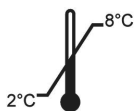

**RUO**



BioGX  
1500 First Avenue, North, L136, Birmingham, AL 35203, USA  
Phone: +1.205.250.8055  
Fax: +1.205.449.8055

## **For Research Use Only**

*Research use only reagents are not intended for human or animal diagnostic use. It is the responsibility of the end user to determine the performance of the reagents in an appropriately designed validation study for their intended use.*

The Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci real-time PCR-based detection reagent is manufactured and packaged as an open system reagent (OSR) for use with open system platforms and has to be validated by the user. Examples of open system platforms are the BD MAX™ System (Windows® Software V4.72A or later), Applied Biosystems QuantStudio™ 5 (Design & Analysis software version 1.5.1 or later), and Bio-Rad CFX96 Touch™ (Maestro software version 1.1 or later) real-time PCR platforms.

**PLEASE READ ENTIRE PACKAGE INSERT BEFORE PROCEEDING TO USE THE OSR.**

### **PRODUCT OVERVIEW**

The BioGX Sample-Ready™ OSR has been formulated in lyophilized format for the multiplex real-time PCR-based detection of DNA from Mycoplasma pneumoniae (Mp181 - CARDS Toxin gene<sup>1</sup>), Legionella spp. (L. pneumophila and L. longbeachae; ssrA gene<sup>1,2</sup>), Chlamydia pneumoniae (argR gene<sup>1</sup>), Chlamydia psittaci (ompA gene<sup>3</sup>) and a Sample Processing Control (SPC) (Drosophila DNA) for use with BD MAX extraction kit. The Drosophila DNA serves as both a sample processing control and an internal amplification control. The configuration compatible with the ABI and Bio-Rad platform targets the human RNase P gene to serve as an endogenous extraction control. The two different formats for the lyophilized Sample-Ready OSR kits are available:

1. BD MAX™ System  
**REF 450-008-C-MAX**
2. ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ Platforms  
**REF 450-008-LMP**

# Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

## **Note:**

**BD MAX™ System OSR (450-008-C-MAX)** contains all PCR primers, probes, enzymes, dNTPs, MgCl<sub>2</sub>, buffers, and other components required for PCR reaction. BD MAX extraction kits available from Becton, Dickinson and Company include a Drosophila DNA sequence which serves as a Sample Processing Control (SPC). No exogenous addition of SPC is needed.

**OSR for ABI and Bio-Rad Platforms (450-008-LMP)** contain all PCR primers, probes, enzymes, dNTPs, MgCl<sub>2</sub>, buffers, and other components required for the PCR reaction. No exogenous addition of SPC is needed.

## **PACKAGE CONTENTS**

**BioGX REF:** 450-008-C-MAX

**Platform(s):** BD MAX™ System

Each 24-reaction package contains two pouches:

1. First pouch contains 24 sealed BD MAX™ tubes of BioGX lyophilized Sample-Ready™ OSR, each tube sufficient for a 12.5 µL PCR reaction.
2. Second pouch contains 24 sealed BD MAX™ tubes, each containing 25 µL of OSR-specific BioGX Rehydration Buffer.

**BioGX REF:** 450-008-LMP

**Platform(s):** ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™

Each 64-reaction package consists of one pouch:

The pouch contains 2 x 8-tube strips. Each tube contains BioGX lyophilized Sample-Ready™ OSR sufficient for 4 x 15 µL PCR reactions.

**Note:** 8-tube strips containing reagents **ARE NOT** compatible with ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ instruments.

## Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

### EQUIPMENT AND MATERIALS REQUIRED BUT NOT PROVIDED

- Lyophilized Positive Control Template DNA Beads ( $10^5$  copies/bead)
  - BioGX Mycoplasma pneumoniae (Part number 720-0010)
  - BioGX Legionella pneumophila (Part number 720-0011)
  - BioGX Legionella longbeachae (Part number 720-0150)
  - BioGX Chlamydia pneumoniae (Part number 720-0012)
  - BioGX Chlamydia psittaci (Part number 720-0149)
- Lyophilized Positive Control Template DNA Beads ( $10^5$  copies/bead) for ABI and Bio-Rad platforms
  - BioGX RNaseP (DNA) (Part number 720-0009)
- BioGX Molecular Grade Water or equivalent
  - BioGX Rehydration Water (Part number: 800-0035-12)
- BioGX BioGX Lyophilized Far Red Dye Calibration Kit for Applied Biosystems QuantStudio 5 (Part number 720-0219)
- BD MAX™ ExK™ DNA-1 (US BD catalog no. 442817/ International BD catalog no. 442818)
- BD MAX™ PCR Cartridges (US and International BD catalog no. 437519)
- Vortex Genie 2 Vortexer (VWR catalog no. 58815-234) or equivalent
- Disposable nitrile gloves
- Applied Biosystems QuantStudio 5 (0.2 mL) consumables
  - Thermo Fisher optical 8-tube strip (catalog no. 4316567)
  - Thermo Fisher ultra-clear optical caps, strips of 8 (catalog no. AB-0866)
  - Thermo Fisher 96-well optical clear reaction plates (catalog no. A36924)
  - Thermo Fisher 96-well qPCR plate seals (catalog no. AB-1170)
- Bio-Rad CFX96 Touch consumables
  - Bio-Rad 8-tube PCR strips without caps (catalog no. TLS0851)
  - Bio-Rad Optical flat 8-cap strips for PCR tubes (catalog no. TCS0803)
  - Bio-Rad 96-well plates (catalog no. HSP9655)
  - Bio-Rad 96-well plate sealing film, optical (catalog no. MSB1001)

### **WARNINGS AND PRECAUTIONS**




- For research use only. Not intended for human or animal diagnostics use.
- If handling biological samples, including used Extraction Kits and PCR Cartridges, treat as if capable of transmitting infectious agents in accordance with safe laboratory procedures such as those described in CLSI Document M29<sup>4</sup> and in Biosafety in Microbiological and Biomedical Laboratories<sup>5</sup>.
- BioGX REF: 450-008-C-MAX has been quality control tested only with the BD Open System Extraction Kits on the BD MAX<sup>TM</sup> System.
- BioGX REF: 450-008-LMP has been quality control tested only with the Applied Biosystems QuantStudio<sup>TM</sup> 5 and Bio-Rad CFX96 Touch<sup>TM</sup> real-time PCR platforms.
- Do not use the reagents if the protective pouches are open or torn upon arrival.
- Close reagent protective pouches promptly with the zip seal after each use. Remove any excess air in the pouches prior to sealing and store at 2-8 °C.
- Do not remove desiccant from the PCR master mix pouches.
- Do not use Sample-Ready<sup>TM</sup> master mix if the desiccant is not present or is broken inside the Sample-Ready<sup>TM</sup> master mix pouches.
- Do not use reagent tubes if the foil seal has been opened or damaged.
- Do not mix reagents from different pouches and/or kits and/or lots.
- Do not use expired reagents and/or materials.
- Refer to BD MAX<sup>TM</sup> ExK<sup>TM</sup> DNA-1 Extraction Kit Instructions or to other respective nucleic acid extraction kits manufacturers' instructions for information about proper handling, cautions, and proper waste disposal.
- Do not mix septum caps between Sample Buffer Tubes or re-use septum caps as contamination may occur and compromise test results.
- Check BD Unitized Reagent Strips for proper liquid fills (ensure that the liquids are at the bottom of the tubes).
- Do not pipette by mouth.
- Do not smoke, drink, or eat in areas where samples or kits are being handled.
- Dispose of unused reagents and waste in accordance with country, federal, provincial, state, and local regulations.
- Use clean gloves when handling extraction kit components and PCR reagents and buffer tubes.



# Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

## STORAGE REQUIREMENTS AND RECOMMENDATIONS

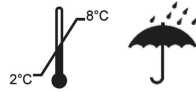

Reagents are stable at a temperature range of 2-30°C during shipment for 5 days, but BioGX recommends long-term storage at 2-8°C. Reagents have been tested to demonstrate optimal performance when stored properly and consumed by the Manufacturer Recommended Use By Date. The end user may opt to extend the useful life for Research Use Only reagents upon completing their own performance validations. BioGX’s guarantee of reagent integrity does not extend beyond the Manufacturer Recommended Use By Date when stored properly. Avoid exposing the reagents (lyophilized or rehydrated) to direct sunlight or long-term ambient lighting. Tightly reseal the pouch with unused reactions and immediately return to a refrigerator after opening. To mitigate reagent performance degradation from exposure to moisture, BioGX suggests using the entire contents of the opened pouch within 1 month; however, the user may choose to verify an extended working time > 1 month by performance testing with positive controls and an examination of the sample preparation control target.

## SAMPLE TYPES UTILIZED FOR QUALITY CONTROL TESTING

- Nasal wash
- Pretreated sputum samples
- Pretreated bronchoalveolar lavage (BAL) samples

## REAGENT OPTICAL CONFIGURATION

**Table 1A.** Optical Channel Configuration for REF 450-008-C-MAX.

Optical Channel (Fluorophore Equivalent)	Target
FAM	Chlamydia psittaci
HEX	Chlamydia pneumoniae
Texas Red	Legionella spp.
Cy5	Mycoplasma pneumoniae
Cy5.5	SPC

## Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

**Table 1B.** Optical Channel Configuration for REF 450-008-LMP<sup>a</sup>.

Optical Channel (Fluorophore Equivalent)	Target
FAM	Chlamydia psittaci
HEX	Chlamydia pneumoniae
Texas Red	Legionella spp.
Cy5	Mycoplasma pneumoniae
Cy5.5 <sup>a</sup>	RNase P

**<sup>a</sup>When using the Applied Biosystems Quantstudio 5 instrument, the BioGX Lyophilized Far Red Dye Calibration Kit for Applied Biosystems QuantStudio 5 (Part number 720-0219) will be needed to add this dye to the instrument dye library. Please refer to the Product Insert of 720-0219 available on [www.biogx.com](http://www.biogx.com) for instructions and further details<sup>2</sup>.**

### **QUALITY CONTROL AND EXTRACTION TESTING PARAMETERS**

#### **BD MAX™ System (BioGX REF: 450-008-C-MAX)**

As a starting point, users can import and install the Electronic User Defined Protocol (eUDP) that utilizes **ExK DNA-1** extraction kits onto the BD MAX™. BioGX's most current eUDP utilized for quality control of this product can be obtained by sending an email to [TS@biogx.com](mailto:TS@biogx.com). Please refer to the BD MAX™ user manual<sup>6</sup> for uploading instructions.

#### **ABI and Bio-Rad Platforms (BioGX REF: 450-008-LMP)**

As a starting point, users can import and install a PCR run file onto:

1. Applied Biosystems QuantStudio™ 5 (Design & Analysis software version 1.5.1 or later)
2. Bio-Rad CFX96 Touch™ (Maestro software version 1.1 or later)

BioGX's most current ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ PCR run files utilized for quality control of this product can be obtained by sending an email to [TS@biogx.com](mailto:TS@biogx.com). Please refer to the Applied Biosystems QuantStudio™5 user manual<sup>7</sup> for uploading instructions. Please refer to the Bio-Rad CFX96 Touch™ user manual<sup>8</sup> for uploading instructions.

---

**SAMPLE VOLUMES UTILIZED DURING QUALITY CONTROL TESTING**

The end user may choose to validate a different pretreatment method or volume of sample to load other than the sample processing used by BioGX for QC testing as outlined below.

**PRETREATMENT PROCEDURES FOR VISCOUS SAMPLES TYPES (i.e., sputum and BAL)**

**APPROPRIATE LOCKING-CAP TUBES OR A LID-LOCK RACK MUST BE USED WHEN SAMPLES ARE BOILED. THE END USER SHOULD USE APPROPRIATE BIOSAFETY PROTOCOLS (INCLUDING A BIOSAFETY HOOD AND RESPIRATOR) WHEN PROCESSING SPUTUM or BAL SAMPLES THAT POTENTIALLY CONTAIN MYCOBACTERIA.**

For viscous samples, the use of a mucolytic agent to reduce viscosity and support efficient DNA extraction is recommended. Examples of three different pretreatment protocols to reduce viscosity of sputum or BAL samples are outlined below. Note: Pretreatment reagents are not included with BioGX 450-008-C-MAX/450-008-LMP.

1. **Proteinase K Pretreatment:** Pellet sample (250 µL sputum or 500 µL BAL) by centrifugation at 20,000 x g for 15 min, decant and wash pellet with 1 mL 20mM Tris HCl pH 8, centrifuge at 20,000 x g for 15 min, decant and add 250 µL of Proteinase K solution (1 mg/mL), incubate at 56°- 65°C for 30 minutes. After Proteinase K digestion, heat to 100°C for 10-15 minutes. Allow for cooling to room temperature.

**BD MAX™ System**

Once the pretreated sample has cooled to ambient temperature, pipette 200 µL of pretreated sample into the BD MAX™ Sample Buffer Tube (SBT), aseptically place the BD™ septum cap on each SBT. Pulse vortex the SBT for 1-3 seconds, and load the SBT into the extraction tray.

**ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ (validated magnetic bead or silica column extraction kits)**

Once the pretreated sample has cooled to ambient temperature, follow manufacturer recommendations for extraction of appropriate sample volume of pretreated sample into validated magnetic-bead or silica column extraction device and proceed with nucleic acid purification.



## **Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents**

2. **Copan SL solution Pretreatment:** (Copan catalog #099CE.A) is a mucolytic agent that supports rapid digestion but does not provide decontamination of natural flora. Manufacturer recommendations should be followed. Depending on the mucopolysaccharide content of the sample, incubation time can range from 15 minutes to 120 minutes.

### **BD MAX™ System**

Pipette 200 µL of pretreated sample into the BD MAX™ Sample Buffer Tube (SBT), aseptically place the BD™ septum cap on each SBT. Pulse vortex the SBT for 1-3 seconds, and load the SBT into the extraction tray.

### **ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ (validated magnetic bead or silica column extraction kits)**

Follow manufacturer recommendations for extraction of appropriate sample volume of pretreated sample into validated magnetic-bead or silica column extraction device and proceed with nucleic acid purification.

3. **BD BBL® MycoPrep™ Pretreatment:** (BD catalog # 240862) supports both mucolytic digestion and decontamination of natural flora. Manufacturer recommendations for BD BBL® MycoPrep™ should be followed. Depending on the mucopolysaccharide content of the sample, incubation time can range from 15 minutes to 30 minutes.

### **BD MAX™ System**

Pipette 200 µL of pretreated sample into the BD MAX™ Sample Buffer Tube (SBT), aseptically place the BD™ septum cap on each SBT. Pulse vortex the SBT for 1-3 seconds, and load the SBT into the extraction tray.

### **ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ (validated magnetic bead or silica column extraction kits)**

Follow manufacturer recommendations for extraction of appropriate sample volume of pretreated sample into validated magnetic-bead or silica column extraction device and proceed with nucleic acid purification.

---

**BioGX QUALITY CONTROL TEST SETUP PROCEDURE**

**BD MAX™ System Quality Control Test Setup**

**Loading a Sample Buffer Tube (SBT)**

1. Add the appropriate sample volume to each SBT.
2. Aseptically place BD™ septum cap on each SBT.
3. Vortex the SBT for 1-3 seconds.
4. Load the SBT into the extraction tray.


**WEAR NITRILE GLOVES WHEN HANDLING LYOPHILIZED REAGENTS TO REDUCE THE GENERATION OF STATIC CHARGES. DO NOT USE LATEX GLOVES.**

**Assembly of BD MAX Extraction Strips and BioGX Reagents**

1. Choose the appropriate BD MAX™ extraction kit (see above). DO NOT use BD MAX™ master mix or the blank 0.3 mL conical tubes that come with the extraction kit.
2. Load the selected extraction cartridges into the extraction tray, 1 per sample to be tested.
3. Snap one BD MAX™ ExK™ DNA Extraction tube into snap-in position 1 (Snap-1) of each extraction strip (**Figure 1**).
4. Snap one BioGX Sample-Ready™ lyophilized reagent tube into position 2 (Snap-2) of each extraction strip. Check to make sure the lyophilized cake is at the bottom of the tube prior to inserting into the strip. The funnel-shaped cake may be in any orientation (v, >, ^, <).
5. Snap one BioGX Rehydration Buffer tube into position 3 (Snap-3) of each extraction strip. Check to make sure the buffer is at the bottom of the tube prior to inserting into the strip. Position 4 (Snap-4) will remain empty.
6. Lift the tray and briefly examine the bottom of each strip to ensure all reagents are at the bottom of each tube and bubbles are not present in Extraction buffer reservoirs.
7. Proceed with worklist generation and sample loading per BD MAX™ operating instructions. Select the appropriate User Defined Protocol (UDP). Load the extraction tray and, if necessary, a new PCR card into the instrument, close the door, and click “Start Run.”
8. Analyze the results by opening the completed run file in the “Results” tab.

# Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

## BD MAX ExK™ 4-Snap Unitized Reagent Strip Single Master Mix Type 4 Setup

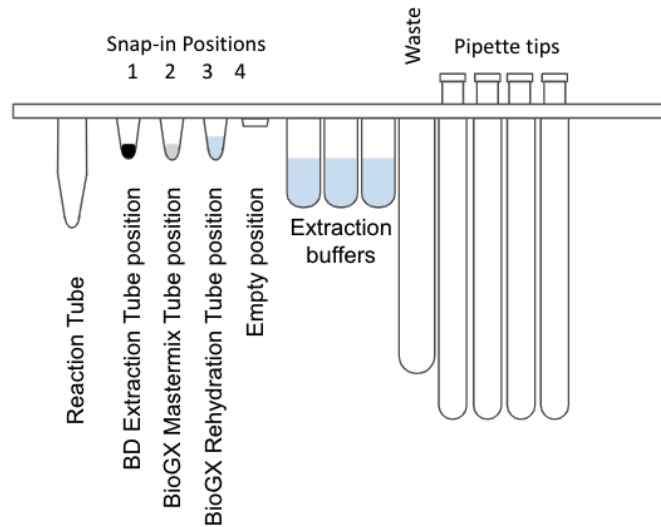

**Figure 1.** Diagram of BD MAX™ ExK™ 4-snap Unitized Reagent Strips.

### Important Note

Always first insert all Snap-1 tubes, then all Snap-2 tubes, then all Snap-3 tubes. The Snap 4 position will remain empty unless the user has set up the reagent to run in dual master mix mode.

Approximately 25 µL of extracted nucleic acid remains in the position 3 tube after extraction. This may be removed and saved for further analyses after the run has been completed.

---

**ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ Quality Control Test Setup**

**WEAR NITRILE GLOVES WHEN HANDLING LYOPHILIZED REAGENTS TO REDUCE THE GENERATION OF STATIC CHARGES. DO NOT USE LATEX GLOVES.**

**Assembly of BioGX Reagents**

1. Prepare the appropriate number of 8-tube PCR strips or 96-well PCR plates.
2. Transfer 40 µL of molecular grade water to one vial of lyophilized BioGX reagents. The rehydrated master mix is more than sufficient for 4 samples to be tested.

**Note:** 8-tube strips containing reagents **ARE NOT** compatible with ABI QuantStudio™ 5 and Bio-Rad CFX96 Touch™ instruments. Rehydrated master mix must be transferred to 8-tube PCR strips or 96-well PCR plates compatible with appropriate instruments.

3. Mix by gently pipetting up and down. (IMPORTANT: Keep rehydrated master mix in a cold block or on ice if setup cannot not be completed within 20 minutes. If the rehydrated master mix cannot be used immediately, it can be capped and stored up to 24 hours at 2-8°C, protected from light).
4. Transfer 10 µL of rehydrated master mix to the bottom of 4 empty wells.
5. To each well containing 10 µL of rehydrated master mix, add 5 µL of extracted sample.
6. Affix the appropriate optical caps or optical plate seals.
7. Pulse spin the sealed PCR plate or tube to mix and bring liquid to the bottom.
8. Load 8-tube PCR strips or 96-well PCR plates into the real-time PCR platform and start the run. Avoid unnecessary delay once tubes/plates are loaded into the real-time PCR instrument.

**ASSAY PERFORMANCE**

All BioGX Research Use Only products are designed to detect 20 copies or less of the target nucleic acid per reaction.

# Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

## INTERPRETATION OF RESULTS

**Table 2A.** Multiplex PCR Results Interpretation for 450-008-C-MAX.

Chlamydia psittaci	Chlamydia pneumoniae	Legionella spp.	Mycoplasma pneumoniae	SPC	Interpretation
+	-	-	-	+/-	Chlamydia psittaci POSITIVE*
-	+	-	-	+/-	Chlamydia pneumoniae POSITIVE
-	-	+	-	+/-	Legionella spp. POSITIVE
-	-	-	+	+/-	Mycoplasma pneumoniae POSITIVE
-	-	-	-	+	Chlamydia psittaci, Chlamydia pneumoniae, Legionella spp., Mycoplasma pneumoniae NEGATIVE
-	-	-	-	-	Unresolved**

\*See *Chlamydia psittaci* clade species and cross-reactivity section below.

\*\*Failed PCR due to inhibition, reagent failure or incorrect assembly of PCR reaction.

**Table 2B.** Multiplex PCR Results Interpretation for 450-008-LMP.

Chlamydia psittaci	Chlamydia pneumoniae	Legionella spp.	Mycoplasma pneumoniae	RNase P	Interpretation
+	-	-	-	+/-	Chlamydia psittaci POSITIVE*
-	+	-	-	+/-	Chlamydia pneumoniae POSITIVE
-	-	+	-	+/-	Legionella spp. POSITIVE
-	-	-	+	+/-	Mycoplasma pneumoniae POSITIVE
-	-	-	-	+	Chlamydia psittaci, Chlamydia pneumoniae, Legionella spp., Mycoplasma pneumoniae NEGATIVE
-	-	-	-	-	Unresolved**

\*See *Chlamydia psittaci* clade species and cross-reactivity section below.

\*\*Failed PCR due to inhibition, reagent failure or incorrect assembly of PCR reaction.

**Chlamydia psittaci clade species and Cross-reactivity**

For samples that are reported in the *Chlamydia psittaci* optical channel, cross-reactivity with the following *Chlamydia* species cannot be ruled out. This assay is not designed for species differentiation of the closely related *Chlamydia psittaci* clade species<sup>10-12</sup> (see list below). Human infection with the below listed organisms is relatively uncommon, but has been reported in the literature<sup>13-15</sup>.

**Chlamydia psittaci clade species**

- *Chlamydia abortus*
- *Chlamydia buteonis*
- *Chlamydia caviae*
- *Chlamydophila felis*

## REFERENCES









1. Thurman, Kathleen A., et al. "Detection of Mycoplasma pneumoniae, Chlamydia pneumoniae, and Legionella spp. in clinical samples using a single-tube multiplex real-time PCR assay." *Diagnostic microbiology and infectious disease* 70.1 (2011): 1-9.
2. Stølhaug, Anne, and Kåre Bergh. "Identification and differentiation of Legionella pneumophila and Legionella spp. with real-time PCR targeting the 16S rRNA gene and species identification by mip sequencing." *Applied and environmental microbiology* 72.9 (2006): 6394-6398.
3. Branley, J. M., et al. "Real-time PCR detection and quantitation of Chlamydia psittaci in human and avian samples from a veterinary clinic cluster." *European Journal of Clinical Microbiology & Infectious Diseases* 27.4 (2008): 269-273.
4. Clinical and Laboratory Standards Institute. Protection of laboratory workers from occupationally acquired infections; Approved Guideline. Document M29 (Refer to the latest edition).
5. Centers for Disease Control and Prevention and National Institutes of Health. Biosafety in Microbiological and Biomedical Laboratories. Choosewood L.C. and Wilson D.E. (eds) (2009). HHS Publication No. (CDC) 21-1112.
6. BD MAX™ System User's Manual (refer to the latest revision) BD Life Sciences, Sparks, Maryland 21152 USA.
7. QuantStudio™ Design and Analysis software User Guide, ThermoFisher Scientific, Waltham, Massachusetts, USA (Refer to the latest version).
8. CFX96 Touch, CFX96 Touch Deep Well, CFX Connect, and CFX384 Touch Systems Instruction Manual (refer to the latest version), Bio-Rad Laboratories, Inc., Hercules, California, USA (Refer to the latest version).
9. BioGX Lyophilized Far Red Dye Calibration Kit for Applied Biosystems QuantStudio™ 5 Real-Time PCR system Product number series 720-0219 ([www.biogx.com](http://www.biogx.com))
10. Longbottom, David, et al. "Whole genome de novo sequencing and comparative genomic analyses suggests that Chlamydia psittaci strain 84/2334 should be reclassified as Chlamydia abortus species." *BMC genomics* 22.1 (2021): 1-18.
11. Everett, Karin DE, Robin M. Bush, and Arthur A. Andersen. "Emended description of the order Chlamydiales, proposal of Parachlamydiaceae fam. nov. and Simkaniaceae fam. nov., each containing one monotypic genus, revised taxonomy of the family Chlamydiaceae, including a new genus and five new species, and standards for the identification of organisms." *International Journal of Systematic and Evolutionary Microbiology* 49.2 (1999): 415-440.
12. Sachse, Konrad, et al. "Emendation of the family Chlamydiaceae: proposal of a single genus, Chlamydia, to include all currently recognized species." *Systematic and applied microbiology* 38.2 (2015): 99-103.
13. Lutz-Wohlgroth, Leslie, et al. "Chlamydiales in guinea-pigs and their zoonotic potential." *Journal of Veterinary Medicine Series A* 53.4 (2006): 185-193.
14. Heddema, Edou R., et al. "Typing of Chlamydia psittaci to monitor epidemiology of psittacosis and aid disease control in the Netherlands, 2008 to 2013." *Eurosurveillance* 20.5 (2015): 21026.
15. Ramakers, Bart P., et al. "Zoonotic Chlamydia caviae presenting as community-acquired pneumonia." *New England Journal of Medicine* 377.10 (2017): 992-994.

**Please call BioGX or email [info@biogx.com](mailto:info@biogx.com) with any questions you may have regarding this product.**

# Mycoplasma pneumoniae, Legionella spp., Chlamydia pneumoniae, Chlamydia psittaci Open System PCR Reagents

Rev. #	Effective Date	Summary of Changes
05	26 JAN 2024	Clarification of species detection in Legionella spp. channel (Legionella pneumophila and Legionella longbeachae). Correction of clerical error.
04	30 SEP 2022	Addition of cross-reactivity statement for C. psittaci primer and probe set
03	25 JAN 2022	Update to include reference to BioGX Far Red Dye Calibration Kit for use with Applied Biosystems Quantstudio 5
02	04 AUG 2021	Update branding, extraction control for ABI and Bio-Rad platforms and storage conditions.
01	28 MAY 2021	Initial Release.

## SYMBOLS

Symbol	Meaning	Symbol	Meaning
	Catalog number		Contains sufficient for <n> tests
	Research Use Only		Manufacturer
	Keep dry		Temperature limitation
	Consult instructions for use		Biological Risks



**BioGX**

BioGX  
1500 First Avenue, North, L136, Birmingham, AL 35203, USA  
Phone: +1.205.250.8055  
Fax: +1.205.449.8055